

firesafe™

safer oxygen therapy



4	EN	Firesafe™ Cannula Valve - Instructions for Use
5	DE	Firesafe™ Rückschlagventil – Bedienungsanleitung
6	FR	Vanne de canule Firesafe™ - Mode d'emploi
7	ES	Válvula para cánulas Firesafe™ - Instrucciones de uso
8	IT	Valvola per cannula Firesafe™ - Istruzioni per l'uso
9	NL	Firesafe™ Cannula Valve - Gebruiksaanwijzing
10	PT	Válvula de cânula Firesafe™ - Instruções de utilização
11	SV	Firesafe™ Slangventil - Användarinstruktioner
12	CZ	Kanylový ventil Firesafe™ - Pokyny k použití
13	PL	Zawór kaniuli Firesafe™ - Instrukcja obsługi
14	RU	Клапан для каниули Firesafe™ - Инструкции по эксплуатации
15	RO	Valva pentru canulă Firesafe™ - Instrucțiuni de utilizare
16	FI	Firesafe™ kanyyliventtiili - Käyttöohje
17	DA	Firesafe™ Cylinderventil - Brugsvejledning
18	SK	Kanylový ventil Firesafe™ - Návod na použitie
19	NO	Firesafe™ kanyleventil - Bruksanvisning
20	SL	Ventil Firesafe™ za kanilo - Navodila za uporabo
21	LT	Firesafe™ kaniulės vožtuvas - Naudojimo instrukcija
22	LV	Firesafe™ kanulas vārsts - Lietošanas instrukcija
23	HU	Használati útmutató Firesafe™ - kanülszelep
24	EL	Βαλβίδα Κάννουλας Firesafe™ - Οδηγίες Χρήσης
25	ET	Firesafe™ kanüüli klapp - Kasutusjuhend
26	BG	Каниопа с клапан Firesafe™ - Инструкции за употреба
27	HR	Ventil kanile Firesafe™ - Upute za uporabu

EN	Figure 1	PT	Figura 1	FI	Kuvan 1	LV	Attēla 1
DE	Abbildung 1	SV	Figur 1	DA	Figur 1	HU	Ábra 1
FR	Figure 1	CZ	Obrázku 1	SK	Obrázku č. 1	EL	Σχήμα 1
ES	Figura 1	PL	Rysunku 1	NO	Figur 1	ET	Joonisel 1
IT	Figura 1	RU	Рисунке 1	SL	Sliki 1	BG	Фигура 1
NL	Afbeelding 1	RO	Figura 1	LT	Pav 1	HR	Slici 1

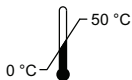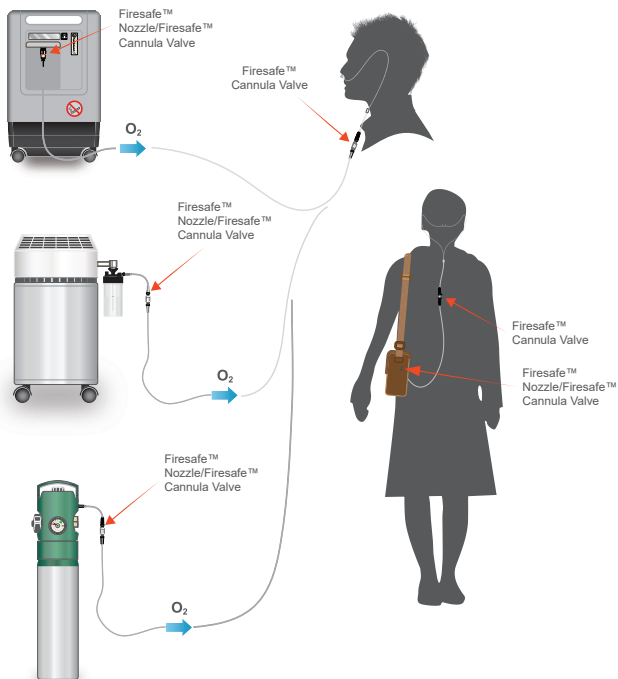EN	Figure 2	PT	Figura 2	FI	Kuvan 2	LV	Attēla 2
DE	Abbildung 2	SV	Figur 2	DA	Figur 2	HU	Ábra 2
FR	Figure 2	CZ	Obrázku 2	SK	Obrázku č. 2	EL	Σχήμα 2
ES	Figura 2	PL	Rysunku 2	NO	Figur 2	ET	Joonisel 2
IT	Figura 2	RU	Рисунке 2	SL	Sliki 2	BG	Фигура 2
NL	Afbeelding 2	RO	Figura 2	LT	Pav 2	HR	Slici 2

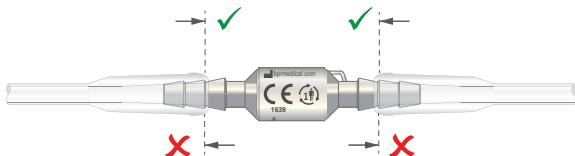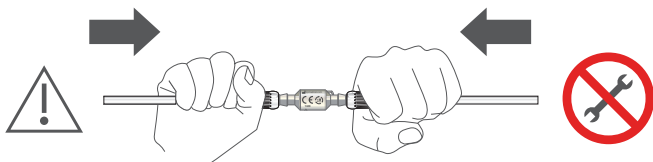Intended Use

The Firesafe™ Cannula Valve is a thermal fuse designed to extinguish an oxygen delivery tube fire by stopping the flow of oxygen if the tube is accidentally ignited and the fire burns back and reaches the device.

Oxygen fires are relatively common, often caused by patients smoking whilst receiving oxygen therapy. Such fires are particularly prevalent where oxygen therapy is provided in the home (domiliary oxygen), where there is limited patient supervision, though they may occur in more supervised environments like healthcare facilities and nursing homes too.

Oxygen is not flammable but oxygen enrichment will dramatically increase the rate and severity of combustion; materials that will not burn in air may do so in an oxygen enriched environment. Should a fire start in the oxygen delivery tube, it will burn back towards the oxygen source, typically an oxygen cylinder, oxygen concentrator or liquid oxygen supply.

It is recommended that two Firesafe™ devices are fitted to each oxygen delivery circuit. The first should be positioned close to the source of supply and the second close to the patient as shown in figure 1. The Firesafe™ Cannula Valve provides no protection if a fire burns through the oxygen delivery tubing upstream of the device.

The Firesafe™ Cannula Valve supplied with this instruction is bidirectional in operation and can be installed either way around in the oxygen delivery tubing.

The Firesafe™ Cannula Valve is intended to be fitted as a risk control measure to low flow, low pressure oxygen delivery circuits for patients requiring supplemental oxygen, whether in the home, healthcare facility or in the care of the emergency services. It may be used for patients for whom supplemental oxygen is critical as long as the potential risks of introducing the device are assessed and additional risk control measures, where necessary, are applied. Risk controls may include monitoring and increased supervision of the patient.

The Firesafe™ Cannula Valve is contraindicated for use downstream of aerosolised drug therapy devices. The device is a non-sterile, single patient use multiple use device.

Warnings!

1. Read through this instruction before installing a Firesafe™ Cannula Valve. As with all medical equipment, attempting to use or install this device without a thorough understanding of its operation and limitations may result in patient or user injury.
2. This device is intended for use in oxygen delivery tubing that has a low static operating pressure not exceeding that stated in the device specification. It must not be used for other applications.
3. This device offers a resistance to flow; consider the impact of this impedance on the system when installing the device.
4. Do not store or install this device near an open flame or near a source of excessive heat that is likely to exceed that stated in the Device Specification.
5. Keep this device free from oil and/or grease. In particular, do not lubricate the hose barb connectors.
6. Never administer oxygen or undertake oxygen therapy while smoking or when near an open flame.
7. Additional risk control measures (e.g. SpO2 monitoring with an alarm) are essential for patients who might suffer injury or death in the event that their oxygen flow is stopped unexpectedly.
8. The device is intended for single patient multiple use and should not be moved between installations due to the risk of cross-contamination.
9. This device is designed for use in the delivery of oxygen or oxygen enriched air. Do not use the device with any other gas.

Cleaning, Maintenance & Disposal

Clean the exterior surfaces of the device using an alcohol or disinfectant wipe. Do not immerse the device in any fluid or allow fluid to enter through the hose barb connections. The Firesafe™ Cannula Valve is maintenance free and has an intended life of 5 years. Dispose of any Firesafe™ Cannula Valve when it reaches the expiry date stated on the device label. Once actuated, the Firesafe™ Cannula Valve cannot be reset and must be discarded.

Device Specification

Resistance to Flow	0.07 kPa at 2 l/min
	0.29 kPa at 5 l/min
Flow Rate ⁽¹⁾	0.97 kPa at 10 l/min
	0.25 - 25 l/min
Maximum Operating Pressure ⁽²⁾	600 kPa
Maximum Internal Leakage Following Activation	10 ml/min
Maximum External Leakage Following Activation	5 ml/min
Operating Temperature	0 to 50 °C (32 to 122 °F)
Transit & Storage Temperature range	-25 to 70 °C (-13 to 158 °F)
Humidity Range	0 to 100% RH

⁽¹⁾The Firesafe™ Cannula Valve can be safely used at lower flow rates but may not meet the internal leakage rate specified when activated in the event of a fire.

⁽²⁾Maximum static pressure delivered by the gas supply source.

Installation Instructions

Warning! When actuated, the Firesafe™ Cannula Valve will stop the gas flow. The pressure upstream of the device will rise until it reaches the maximum output pressure of the supply system. In order to ensure that the connection between the Firesafe™ Cannula Valve and the supply system can withstand this pressure, it is recommended that a type test be carried out with each style of connector or tubing used.

Note: Due to the differences in the flexibility of the materials used for tubing and trumpet connectors and variations of inner bore diameter, it may not always be possible to fully engage the second barb as shown in figure 2.

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